


## HydroCheck-E

## Easily test for residual water in flexible endoscopes channels





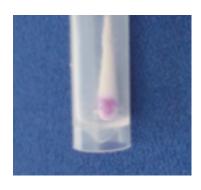

## Product Code:77080-06A (1.7mm), 77080-06B (2.8mm), 77080-06C (3.8mm)

Following the cleaning process great care must be taken to ensure that flexible endoscopes are thoroughly dry before use as the presence of any residual moisture poses the very real risk of micro-organism development within the scope.

The Hydrocheck-E is a simple to use, frictionless, pre-prepared swab test for checking to ensure that flexible endoscope channels are truly dry.

The tests are available with swabs in three sizes; 1.7mm, 2.8mm & 3.8mm

The test is conducted by passing the highly absorbent cotton swab through the channel and then placing the swab into the Hydrocheck-E test vial.

The presence of the smallest amount of moisture on the swab will produce a purple dot on the swab surface and if large amounts of moisture are detected the entire vial solution can change colour.

With a detection limit as low as  $0.05 \,\mu$ l & producing immediate test results, the Hydrocheck-E provides confidence in your flexible endoscope re-processing regime.



Web: www.minervastericlean.co.uk Email: info@minervastericlean.co.uk Tel: +44 (0) 1249 820178